

**Brookstone®**

---

**MAX®2 Cordless**  
**POWERFUL PERCUSSION MASSAGE**



## TABLE OF CONTENTS

Warnings and Cautions .....	2
Location of Parts and Controls .....	6
Massager Operation .....	7
Product Care and Maintenance .....	10
Specifications .....	11
Warranty .....	13

## CAUTION

- TO REDUCE THE RISK OF ELECTRIC SHOCK, DO NOT REMOVE COVER. THERE ARE NO SERVICEABLE PARTS INSIDE.
- TO REDUCE THE RISK OF FIRE OR ELECTRIC SHOCK, DO NOT EXPOSE THIS UNIT TO RAIN OR MOISTURE.



The lightning flash with arrow-head symbol within an equilateral triangle is intended to alert the user to the presence of uninsulated “dangerous voltage” within the unit’s enclosure that may be of sufficient magnitude to constitute a risk of electric shock.



The exclamation point within an equilateral triangle is intended to alert the user to the presence of important operating and maintenance (servicing) instructions in the literature accompanying the unit.

## IMPORTANT SAFETY INSTRUCTIONS

All of the safety and operating instructions should be read, adhered to and followed before the unit is operated.

## SAVE THESE INSTRUCTIONS

## **DANGER**

To reduce the risk of electric shock, burns, fire or injury:

1. Always unplug the product from the electrical outlet immediately after using and before cleaning.
2. Do not use while bathing or in a shower.
3. Do not place or store product where it can fall or be pulled into a tub or sink.
4. Do not place in or drop into water or other liquid.
5. Do not reach for product that has fallen into water. Unplug it immediately.

## **WARNING**

To reduce the risk of burns, fire, electric shock or injury to persons:

1. Do not leave product unattended when plugged in. Turn off and unplug from outlet when not in use.
2. Close supervision is necessary when this appliance/product is used by or near children or mentally disabled individuals.
3. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water.

4. Keep cord away from heated surfaces.
5. Never block the air openings of the product or place it on a surface, such as a bed or couch, where the air openings may be blocked. Keep the air openings free of lint and hair.
6. Never use while sleeping.
7. Never drop or insert objects into any openings.
8. Do not operate where aerosol products are being used or where oxygen is being administered.
9. Do not use on or near eyes or other highly sensitive areas.
10. Do not operate under a blanket or pillow. Excessive heating can occur and cause fire, electric shock or injury.
11. Do not carry this product by AC Adapter cord or use cord as a handle.
12. Do not use directly on one concentrated area of the body. It is meant to be used in a continuous motion, constantly changing the massage area.
13. To disconnect, turn all controls to the off position and then remove the plug from the outlet.
14. Do not perform any servicing on the product.
15. This product is designed for household use only. Do not use outdoors.
16. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.

## CAUTION



1. Consult your doctor before using this product.
2. Individuals with pacemakers should not use this product.
3. Never leave the product operating or unattended, especially if children are present.
4. Never use the product on open skin eruptions, swollen or inflamed areas.
5. Do not operate after taking pain killers, sedatives or alcohol. Consult your doctor if you are under these medications.
6. Do not operate for more than 20 minutes continuously. It is recommended to turn off the product and allow it to rest (cool) for 40 minutes before re-using. This will prolong the life of the product.
7. The product is designed as a non-professional massager to soothe tired and aching muscles. It is not intended to be a substitute for medical attention. If your ailments and symptoms persist or worsen after use, contact your doctor immediately.

**FOR ADULT USE ONLY.** Remove all scarves, neckties, necklaces, and jewelry before using this massager.

**IMPORTANT:** Any individual who may be pregnant, has a pacemaker, suffers from diabetes, phlebitis and/or thrombosis, is at an increased risk of developing blood clots, or who has pins/screws/artificial joints or other medical devices implanted in his/her body should consult with a physician before using a massaging device designed for home use.

## LOCATION OF PARTS AND CONTROLS

1. Adapter Jack
2. Handle 1
3. On/Off Button
4. Speed Button
5. Speed/Charging LEDs
6. Programs Button
7. Programs LEDs
8. Handle 2
9. Massaging Nodes (2)
10. AC Adapter

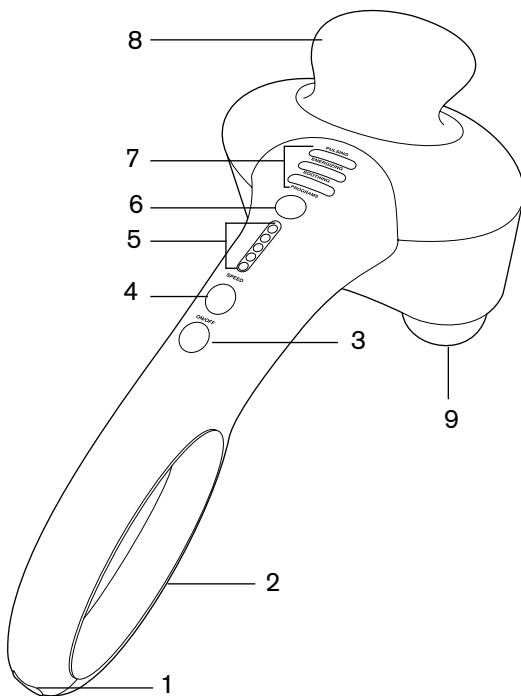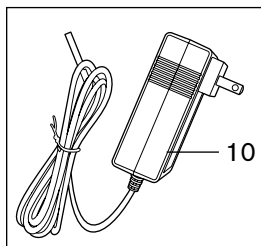

# MASSAGER OPERATION

## CHARGING THE MASSAGER

**Important!** Fully charge the massager for 2 hours before first use.

1. To charge the massager, connect the ADAPTER to the massager and plug it into a wall outlet.
2. The one through five SPEED LEDs will illuminate sequentially while charging, and they will all remain illuminated when charging is complete.
3. A completely depleted massager battery will take up to 2 hours to charge.
4. The average run-time for the battery at the highest setting is 3 hours.
5. The SPEED LEDs will flash green when the massager's battery needs to be recharged.

**Note:** You may also use the massager while it is plugged in to the electrical outlet, however, the massager's battery will not charge while in use.

## TURNING POWER ON AND OFF

1. Press POWER. When the massager is first turned on, all five green SPEED LED lamps light up, from Speed 1 to Speed 5. The blue Program lamps (SOOTHING, ENERGIZING and PULSING) light up once sequentially. The motor then runs at its top speed (Speed 5).
2. Press POWER again to shut off. The motor stops and all SPEED LED lamps shut off.

**Note:** The massager will automatically shut off after 20 minutes.



## USING SPEED

Because the massager runs at Speed 5 when it is first powered on, SPEED must be pressed to change speeds.

- Press SPEED once and the motor turns to the minimum, Speed 1. One green speed-indicating LED lamp lights up.
- Press SPEED again and the motor turns to Speed 2. Speed 1 and Speed 2 green LED lamps light up.
- Press SPEED again and the motor turns to Speed 3. Speed 1, Speed 2 and Speed 3 green LED lamps light up.
- Press SPEED again and the motor turns to Speed 4. Speed 1, Speed 2, Speed 3 and Speed 4 green LED lamps light up.
- Press SPEED again and the motor turns to Speed 5. Speed 1, Speed 2, Speed 3, Speed 4 and Speed 5 green LED lamps light up.
- Press SPEED again and the motor returns to Speed 1.

**Note:** Pressing SPEED thereafter will continuously change the speed sequentially, as described above.

## USING PROGRAMS

- Press the PROGRAMS button once and PULSING is activated.  
The corresponding blue LED lamp lights up. The motor continuously runs in sequence of Speed 5 and a brief stop. The corresponding Speed 5 green lamp lights up.
- Press the PROGRAMS button again and ENERGIZING is activated.  
The corresponding blue LED lamp lights up. The motor continuously runs in sequence of Speed 1, Speed 2, Speed 3, Speed 4; Speed 1, Speed 2, Speed 3, Speed 4 and so on. The corresponding green Speed lamps light up.
- Press the PROGRAMS button again and SOOTHING is activated.  
The corresponding blue LED lamp lights up. The motor continuously runs in sequence of Speed 1, Speed 2, Speed 3, Speed 1, Speed 2 and so on.  
The corresponding green Speed lamps light up.

### Note:

- Pressing PROGRAMS thereafter will continuously change the programs sequentially, as described above.
- Whenever PROGRAMS is exited by pressing SPEED, the massager returns to Speed 5 as the default.

## **PRODUCT CARE AND MAINTENANCE**

1. Place the massager in a safe, cool and dry place when not in use.
2. Store the Adapter with the massager.
3. Wipe massager with a soft, damp cloth using no abrasive cleaner.
4. Never submerge the unit in any liquid.
5. Keep away from all solvents and harsh detergents.
6. DO NOT attempt to repair this massager. There are no user-serviceable parts.
7. Dispose of this product as you would a battery. Contact your local recycling center for disposal rules.

# SPECIFICATIONS

## Adapter

Input	100-240 ~ 50/60Hz 1.5A
Output	15VDC, 3.0A

## RPM

Speed Setting	RPM	RPM Tolerance
1	1,750	+/- 10%
2	2,050	+/- 10%
3	2,250	+/- 10%
4	2,450	+/- 10%
5	2,550	+/- 10%

# NOTES

# ONE (1) YEAR LIMITED WARRANTY

Brookstone warrants this product against defects in materials and/or workmanship under normal use for a period of ONE (1) YEAR from the date of purchase by the original purchaser ("Warranty Period"). If a defect arises and a valid claim is received within the Warranty Period, at its option, Brookstone will either 1) repair the defect at no charge, using new or refurbished replacement parts, or 2) replace the product with a new product that is at least functionally equivalent to the original product, or 3) provide a store credit in the amount of the purchase price of the original product. A replacement product or part, including a user-installable part installed in accordance with instructions provided by Brookstone, assumes the remaining warranty of the original product. When a product or part is exchanged, any replacement item becomes your property and the replaced item becomes Brookstone's property. When a store credit is given, the original product must be returned to Brookstone and becomes Brookstone's property.

**Obtaining Service:** To obtain warranty service, call Brookstone Limited Warranty Service at 1-800-292-9819. Please be prepared to describe the product that needs service and the nature of the problem. A purchase receipt is required. All repairs and replacements must be authorized in advance. Service options, parts availability and response times will vary. You are responsible for delivery and the cost of delivery of the product or any parts to the authorized service center for replacement, per our instructions.

**Limits and Exclusions:** Coverage under this Limited Warranty is limited to the United States of America, including the District of Columbia and the U.S. Territories of Guam, Puerto Rico, and the U.S. Virgin Islands. This Limited Warranty applies only to products manufactured for Brookstone that can be identified by the "Brookstone" trademark, trade name, or logo affixed to them or their packaging. The Limited Warranty does not apply to any non-Brookstone products. Manufacturers or

suppliers other than Brookstone may provide their own warranties to the purchaser, but Brookstone, in so far as permitted by law, provides these products "as is." This warranty does not apply to: a) damage caused by failure to follow instructions relating to product's use or the installation of components; b) damage caused by accident, abuse, misuse, fire, floods, earthquake or other external causes; c) damage caused by service performed by anyone who is not a representative of Brookstone; d) accessories used in conjunction with a covered product; e) a product or part that has been modified to alter functionality or capability; f) items intended to be periodically replaced by the purchaser during the normal life of the product including, without limitation, batteries or light bulbs; g) any product sold "as is" including, without limitation, floor demonstration models and refurbished items; or h) a product that is used commercially or for a commercial purpose.

**BROOKSTONE SHALL NOT BE LIABLE FOR INCIDENTAL OR CONSEQUENTIAL DAMAGES RESULTING FROM THE USE OF THIS PRODUCT, OR ARISING OUT OF ANY BREACH OF THIS WARRANTY. TO THE EXTENT PERMITTED BY APPLICABLE LAW, BROOKSTONE DISCLAIMS ANY AND ALL STATUTORY OR IMPLIED WARRANTIES, INCLUDING, WITHOUT LIMITATION, WARRANTIES OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE AND WARRANTIES AGAINST HIDDEN OR LATENT DEFECTS. IF BROOKSTONE CANNOT LAWFULLY DISCLAIM STATUTORY OR IMPLIED WARRANTIES, THEN TO THE EXTENT PERMITTED BY LAW, ALL SUCH WARRANTIES SHALL BE LIMITED IN DURATION TO THE DURATION OF THIS EXPRESS WARRANTY.**

Some states disallow the exclusion or limitation of incidental or consequential damages or how long an implied warranty lasts, so the above exclusions or limitations may not apply to you. This warranty gives you specific legal rights and you may also have other rights, which vary from state to state.

Find thousands more  
great ideas online

**Brookstone**<sup>®</sup>

792451

Merrimack, New Hampshire USA 03054 • 800-846-3000 • Brookstone.com

Download from [www.Somanuals.com](http://www.Somanuals.com). All Manuals Search And Download.



## Free Manuals Download Website

<http://myh66.com>

<http://usermanuals.us>

<http://www.somanuals.com>

<http://www.4manuals.cc>

<http://www.manual-lib.com>

<http://www.404manual.com>

<http://www.luxmanual.com>

<http://aubethermostatmanual.com>

Golf course search by state

<http://golfingnear.com>

Email search by domain

<http://emailbydomain.com>

Auto manuals search

<http://auto.somanuals.com>

TV manuals search

<http://tv.somanuals.com>